

CHEMGLASS

# Air-Operated Drain Valves

for 10 through 50L Chemglass Reactors

- Fits existing Chemglass reactors with 2" beaded pipe drains
- Ideal for solvent extractions & partitions after a reaction
- Great for low temp reactions where standard valves need constant tightening to prevent leakage
- Held closed by spring force & opened by air pressure
- Chemglass Part Nos. **CG-1968-V-series**



### Control Panel


Included are two, factory-preset regulators. The single control knob quickly opens and closes the valve as needed. Set-up is easy and flexible as the panel may be mounted on the stand or outside a fume hood.

**Please Note:** Due to the automatic nature of this valve, we do not recommend it for slurries or viscous materials.



### Bottom Valve

Features a 3/4" compression fitting on the drain outlet. Nose and housing are PTFE and the shaft is available in 316L stainless or Hastelloy®. The valve mounts easily on existing 2" beaded pipe drain.

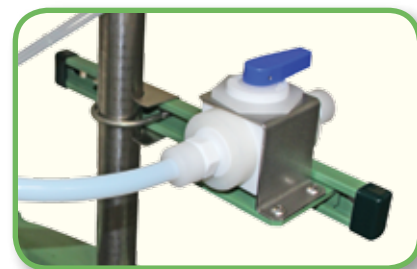

### Optional Secondary Valve Kit

For emergency shut-off of the reactor contents in the unlikely event of valve failure.



See back page  
for ordering  
information. ►

## Air-Operated Drain Valves for 10 through 50L Chemglass Reactors




### Drain Valves with Stainless Steel or Hasteloy® Shafts

Part Number	Description	Price
CG-1968-V-01	Air-Operated Lower Drain Valve, <b>Stainless Steel Shaft</b> , Complete with Air Control Panel, Perfluoro O-rings, Tubing and Fittings	<b>*\$2,950.00</b>
CG-1968-V-02	Air-Operated Lower Drain Valve, <b>Hasteloy® Shaft</b> , Complete with Air Control Panel, Perfluoro O-rings, Tubing and Fittings	<b>*\$3,950.00</b>

Mounting clamp not included. See **CG-9253** Universal Kwik Clamp.

**Please Note:** Due to the automatic nature of this valve, we do not recommend it for slurries or viscous materials.

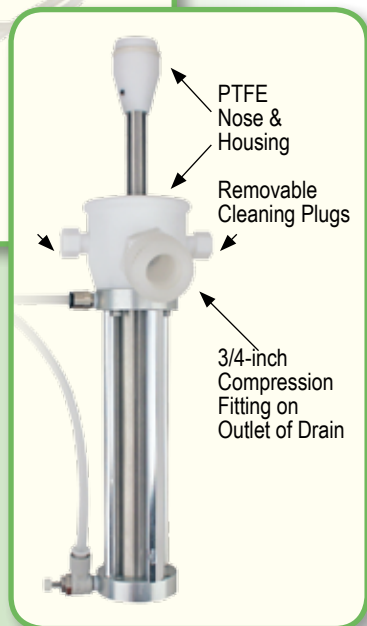

See it in action on VIMEO.com

Search for: Chemglass

It's also on the product page of our website.

Search for: CG-1968-V



### Optional Secondary Valve Kit

Part Number	Description	Price
CG-1968-V-10	Secondary Valve Kit, All-PTFE Ball Valve with 3/4" Compression Fittings for Use with Air-Operated Drain Valves	<b>*\$1,410.00</b>

