

Sterile Sample Tube Assemblies 2026

Custom
Assemblies
available too!

For Sampling and Transfer in Closed Systems

- **over-molded dip tubes** eliminate risk of detachment
- polypropylene centrifuge tubes
- completely **closed cap systems**
- include tubing, pinch clamps and vent filters

Sterile Sample Tube Assemblies with C-Flex® Inlets

Polypropylene centrifuge tube assemblies feature closed cap assemblies with 18" of Size 16 (1/8" x 1/4") C-Flex tubing inlet line with pinch clamps and plugs. They include 3" vent lines and PTFE filters. Certified sterile.

Assemblies can be used for sampling and transfer of biopharmaceutical products and reagents in closed systems. Individually bagged in an outer bag of ten.

Sample Tube Assemblies with C-Flex Dip Tubes, Sterile

Part No.	Vol. (mL)	Cent. Tube Material	Inner Dip Tube Material	Qty.	Price	
CLS-1407-T-312	50	PP	C-Flex	10	\$513.65	*
CLS-1407-T-312A	50	PP	C-Flex	10	\$514.80	*
CLS-1407-T-325	25	PP	C-Flex	10	\$427.25	*


Sterile Sample Tube Assemblies with C-Flex to Tygon® to C-Flex Inlets

Polypropylene centrifuge tube assemblies feature closed cap assemblies with 3/4" of Size 16 (1/8" x 1/4") C-Flex tubing inlet line which transitions to 2 1/2" of Size 14 (1/16" x 1/4") Tygon and then transitions to 17 1/4" of Size 16 (1/8" x 1/4") C-Flex tubing with pinch clamps and plugs. They also include 3" vent lines and PTFE filters. Assemblies are certified sterile.


Assemblies can be used for sampling and transfer of biopharmaceutical products and reagents in closed systems. Individually bagged in an outer bag of ten.

Sample Tube Assemblies with C-Flex Dip Tubes, Sterile

Item No.	Vol. (mL)	Inlet Length (in)	Vent Line Length (in)	Qty.	Price	
CLS-1407-T-310	50	18	3	10	\$570.95	*

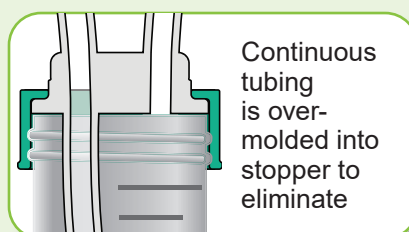

*Life Science Equipment Product Line

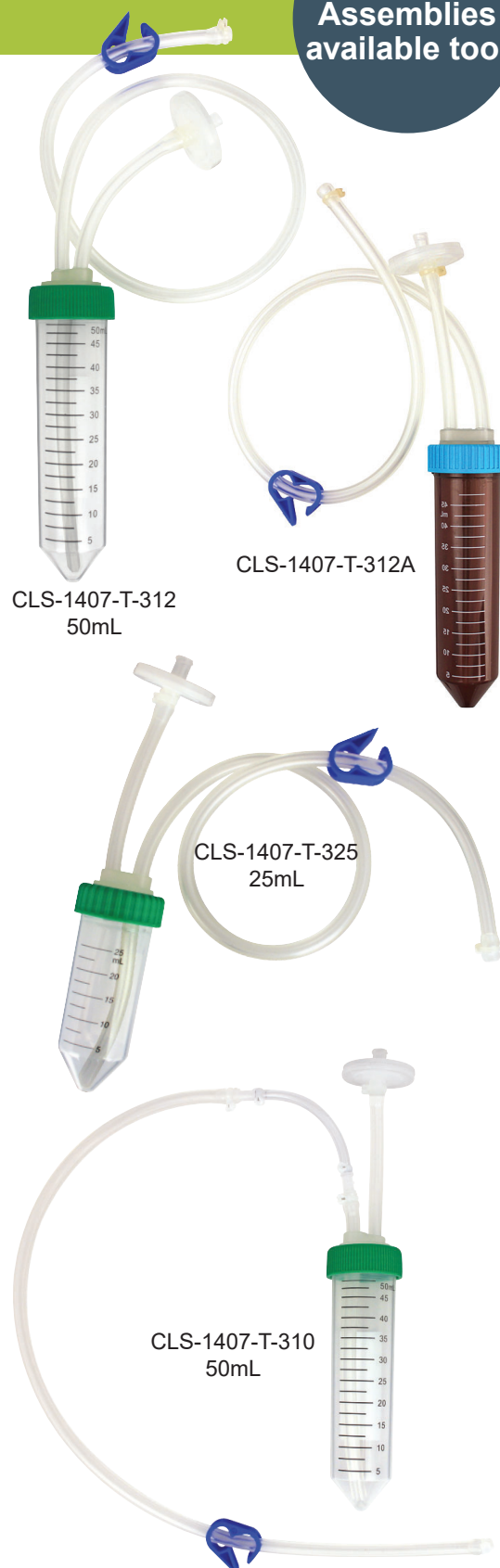

Need help with a custom sample tube size or configuration? Just get in touch.



www.cglifesciences.com - singleuse@cglifesciences.com
T 1.800.843.1794 - F 1.800.922.4361

Sterile Sample Tube Assemblies



CLS-1407-T-515



CLS-1407-T-516

Sample Tube Assemblies, 15mL Centrifuge Tube With TPE and Male Luer Lock Inlet, Sterile

The 15mL polypropylene centrifuge tube assembly utilizes a closed cap assembly with 18" of 1/16" ID x 1/8" OD TPE tubing inlet line with a pinch clamp and a capped male luer lock end. The assembly includes a 3" vent line with a 25mm PES filter. Assembly is certified sterile.

The assembly can be used for sampling and transfer of biopharmaceutical products and reagents in closed systems. Sold individually wrapped in a bag of five or ten

Sample Tube Assemblies with TPE, Male Luer Lock Inlet, Sterile

Part No.	Vol. (mL)	Cent. Tube Material	Inner Dip Tube Material	Qty.	Price
CLS-1407-T-515	15	PP	TPE	5	\$385.75 *
CLS-1407-T-1015	15	PP	TPE	10	\$680.85 *

Sample Tube™ Assemblies, 15mL, Sterile, Single-Use, 15mL Centrifuge Tube

The 15mL polypropylene centrifuge tube assembly utilizes a closed cap assembly with 2" of 1/16" ID x 1/8" OD TPE to 16" of Size 16 (1/8" ID x 1/4" OD) C-Flex tubing inlet line with a pinch clamp and a plugged end. The assembly includes a 3" vent line with a 25mm PES filter. Assembly is certified sterile.

The assembly can be used for sampling and transfer of biopharmaceutical products and reagents in closed systems. Sold individually wrapped in a bag of five or ten

Sample Tube Assemblies with C-Flex Dip Tubes, Sterile

Part No.	Vol. (mL)	Cent. Tube Material	Inner Dip Tube Material	Qty.	Price
CLS-1407-T-516	15	PP	C-Flex	5	\$385.75 *
CLS-1407-T-1016	15	PP	C-Flex	10	\$680.85 *

*Life Science Equipment Product Line

Need help with a custom sample tube size or configuration? Just get in touch.



www.cglifesciences.com - singleuse@cglifesciences.com
T 1.800.843.1794 - F 1.800.922.4361