



Digital Hotplate Stirrer Kits with Tri-Blocks™ and Trial Packs of Vials

Get more than \$200 in **COST SAVINGS!**

Use Promo Code COMBO200



Illuminated icons for Heat, Probe, Setpoint, Timer and Autoreverse



Maintains set RPM even with viscosity changes



5-YEAR WARRANTY



Temperature probe included

Laser-Verified High Performance Control

- Fast, precise temperature control from 0 to 370°C
- Extremely accurate stirring from 30 to 1700rpm
- Auto-reverse option for enhanced mixing
- Torque compensation maintains set RPM as viscosity changes
- Chemical- and scratch-resistant, CerAlTop™, ceramic-coated top

Designed with Advanced Features

- Integral timer - settable from 1 minute to 99 hours, 59 minutes
- Service life meter - track service interval or service lifetime
- Bright white LED display for high-contrast easy reading
- 135mm diameter, ceramic-coated, aluminum alloy plate

Engineered for Safety


- Overheat protection & adjustable safety temperature circuit
- Timer automatically turns heater off after the elapsed time

Watch Video Online



See a case study and more information on the back.

Tri-Block Holds Three Popular Vial Sizes Simultaneously

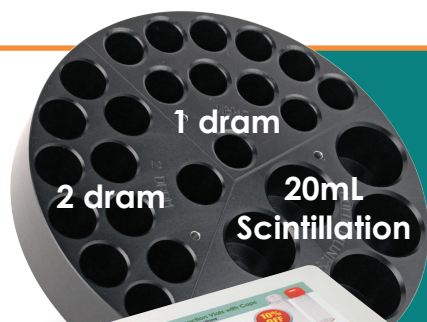

FREE Trial Pack of Pressure Relief Reaction Vials



Tri-Blocks


- Holes spaced for optimal stirring
- Aluminum for efficient heat transfer
- Thermowells provided for temperature probe

Pressure Relief Reaction Vials

- Heat- and chemical-resistant borosilicate glass reaction vials
- Caps prevent pressure from exceeding 150psig
- Trial pack includes 1 and 2 dram vials and 20mL scintillation vials

1-800-843-1794 • www.cglifesciences.com • tech@cglifesciences.com

Fast, Precise Heating and Temperature Control

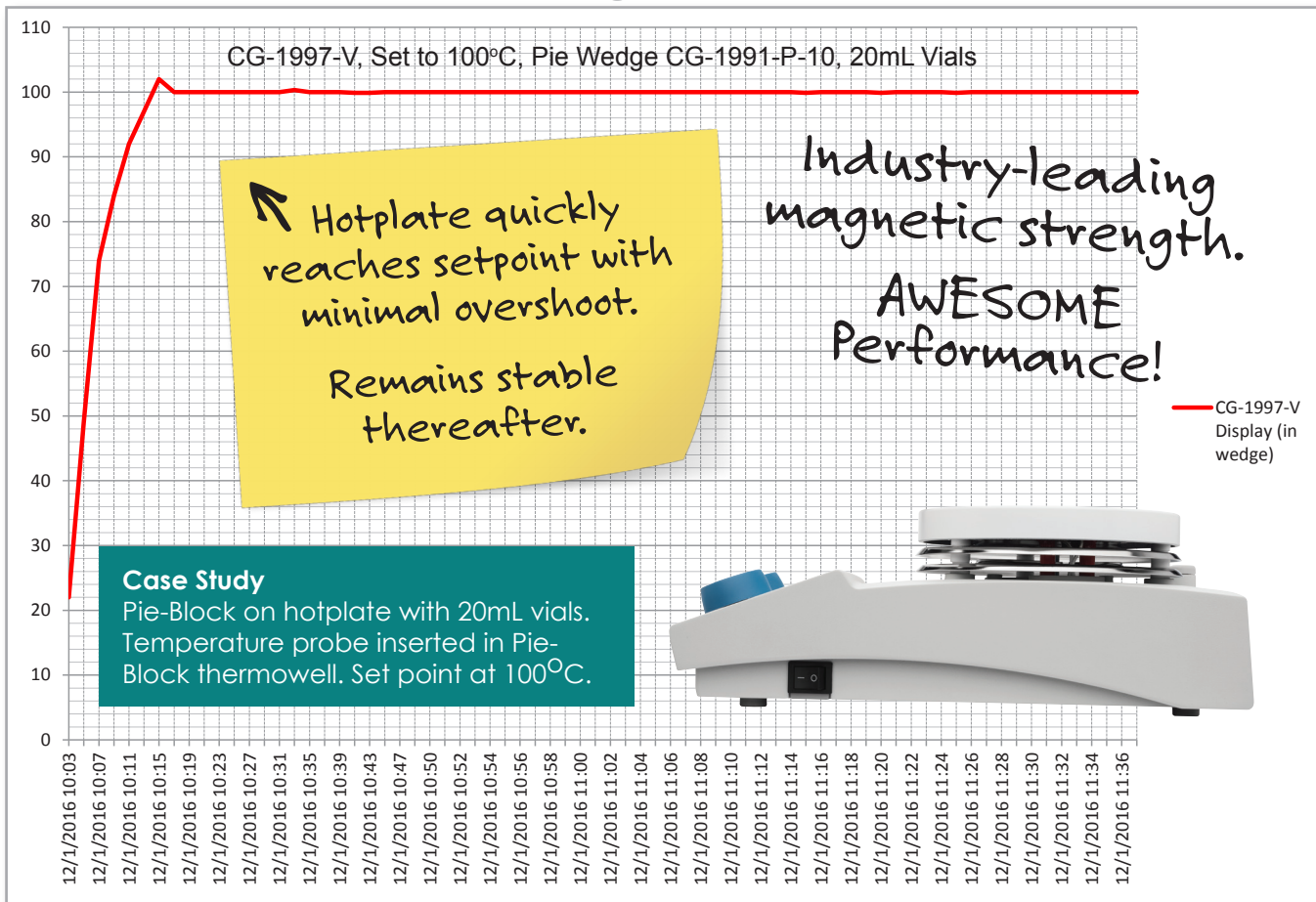

Printed in US. CG-17054R19

CG-1997-V-210

Reg. ~~\$1050~~
\$850*

Tri-Block & Vial Trial Pack Included

CG-1997-V-250

Reg. ~~\$1100~~
\$905*

Tri-Block & Vial Trial Pack Included

Part No.	Description	Price
CG-1997-V-210	Digital Hotplate Stirrer and Tri-Block, AREX-6, Temperature Probe, 115V	\$850.00 *
CG-1997-V-250	Digital Hotplate Stirrer, Tri-Block and Hardware, AREX-6, Temperature Probe, Vertical Support Rod, Boss Head and Support / Holding Rod, 115V	\$905.00 *

1-800-843-1794 • www.cglifesciences.com • tech@cglifesciences.com

RULES - Promo Code COMBO200

Participation in this offer is restricted to purchases at list price of CG-1997-V-210 (COMBO Price of \$825.00) and CG-1997-V-250 (COMBO Price of \$879.00) only. Offer valid from 1/3/2019 through 12/31/2019. FREE Reaction Vial Trial Pack only available with purchase of stirring hotplate and block. Please include promo code COMBO200 when placing the order.

Valid on US and Canadian purchases only. May not be combined with any other offer or discount. This offer is non-transferable and void where prohibited by law or company policy. Offer is subject to Chemglass Life Sciences' standard terms, conditions and contracts. Offer not valid through distributors.

*Not subject to discount.